
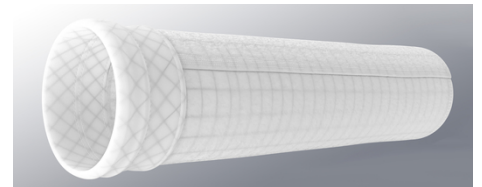


Nutsche Filters

Sefar provides type-specific filter discs and bags for nutsche filters that meet the high standards of the pharmaceutical and chemical industry with regard to certification and traceability.

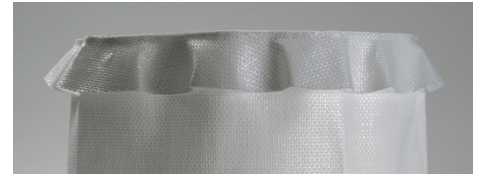

Product Features

Sefar filter solutions for Nutsche Filter

Pressure Nutsche Filters are simple and proven filtration systems for solid-liquid separation of medium-sized batches on various filter media.

In the pharmaceutical and chemical industries, pressure filters and single-layer filters are used in a large variety of forms.

DOWNLOADS

Brochure: Process Filtration (PDF 5300 kb)

Fabrication solutions

Sefar's filter discs are characterized by their precise fit and easy installation. The disc edges can be fabricated using laser cutting, hemmed or taped with reinforcement. In our wide range of filter materials (polypropylene, polyester, PEEK and PTFE fabric), we can supply you with the optimal solution for your application's chemical and thermal requirements. Special bag versions are available as customized solutions.

Unique marking



Customized bag version with flange

Sefar Inc.


Kansas City, MO – USA

Phone 800 995 0531
 Fax +1 816 452 2183

sales.us@sefar.com

[Go to product page](#)

Hemmed with cord

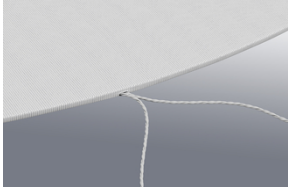

Sefar filter solutions for OEM

Our products meet the specific needs of these filtration applications and are successfully running all over the world on all known nutsche filter brands and OEMs, such as:

Comber Heinkel Mavag AG Waldner

Bachiller Delta

Filter media technology

Sefar's **Double Layer Weave** technology combines fine filter fabric together with a rugged backing cloth.

SEFAR **TETEX** DLW fabrics are produced in our European facilities to meet the highest quality requirements and offer these decisive advantages:

Sefar Inc.

Kansas City, MO – USA

Phone 800 995 0531
Fax +1 816 452 2183

sales.us@sefar.com

[Go to product page](#)

- Low downtime
- High product yield
- Ease of installation
- Excellent cake release
- Minimal product contamination

In addition to DLW fabrics, monofilament and multifilament fabrics are also available to ensure an optimal design.

Double-layer fabric



- With integrated backing cloth for excellent drainage and stability

Monofilament fabric



- For maximum permeability and optimal cleaning properties

Multifilament fabric



- For maximum particle retention, strength and cloth flexibility

Specially made for life science industries

Sefar controls every crucial production step from yarn production to the ready-to-use filters making possible a complete traceability of all critical parts.

Our unique market position in **life science** enables us to supply products to our customers that meet all required quality levels and compliances.

Sefar has perfected its manufacturing processes to meet the needs of the pharmaceutical industry. Our exclusive GMP compliant filter products are manufactured in our cleanroom following strict GMP guidelines.

Sefar Inc.

Kansas City, MO – USA

Phone 800 995 0531
 Fax +1 816 452 2183

sales.us@sefar.com

[Go to product page](#)



- EC 1935/2004
- EU 10/2011
- GMP EC 2023/2006 or cGMP
- FDA CFR 177
- ISO 9001
- Single packaging
- Customized marking
- Traceability



CONTACT

■ Local contact

Sefar Inc.

Kansas City, MO – USA

Phone 800 995 0531
Fax +1 816 452 2183

sales.us@sefar.com

[Go to product page](#)



- [Contact form](#)
- [Send mail](#)

Please call us for further information:
Phone CH: +41 71 898 5700

Locations



Sefar BDH Inc.

200 rue Clement Gilbert
Saguenay QC, G7H 5B1 – Canada
Phone: +1 418 690 0888
Fax: +1 418 690 9499

■ [E-Mail](#)



Sefar Inc.

Kansas City, MO – USA
Phone: 800 995 0531
Fax: +1 816 452 2183

■ [E-Mail](#)



Sefar Inc.

Los Angeles, CA – USA
Phone: 800 995 0531
Fax: +1 909 544 5901

■ [E-Mail](#)

Sefar Inc.

Kansas City, MO – USA

Phone 800 995 0531
Fax +1 816 452 2183

sales.us@sefar.com

[Go to product page](#)



Sefar Inc.

Gray, GA – USA
Phone: 800 995 0531
Fax: +1 478 986 6953

■ [E-Mail](#)



Sefar Inc.

111 Calumet Street
Buffalo, NY 14043 – USA
Phone: 800 995 0531
Fax:

■ [E-Mail](#)

Sefar Inc.

Kansas City, MO – USA

Phone 800 995 0531
Fax +1 816 452 2183

sales.us@sefar.com

[Go to product page](#)